


Oxygen Switch-off Valve with Auxiliary Instant Supply Valve

In critical care units where resuscitation activities require usage of Ambu with instant oxygen supply, SILBERMANN's switch-off valve provides the ideal solution.

The oxygen switch shut off valve is used with an AMBU bag connected to an auxiliary valve. When the bag is lifted from the lever arm the Instant flow valve provides oxygen at a flow rate of up to 25 liter/min. When bag is returned to arm the flow stops immediately. The flowmeter is independent and can be operated regardless the position of the lever arm.

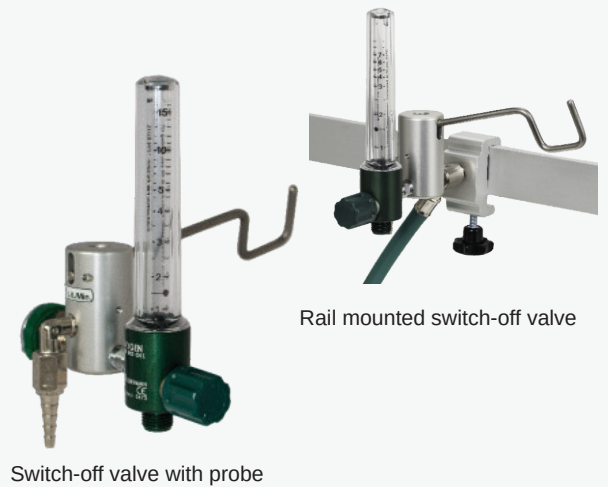

Switch-off valve with probe

Rail mounted switch-off valve



Features & Benefits

- Accurate & reliable gas flow
- Cleaned for oxygen use
- Flowmeter operates independently regardless the position of the lever arm
- Various flow rates available through the auxiliary instant supply valve
- Flowmeter flow rates:
 - Standard: 0-15LPM.
 - Pediatric: 0-7LPM.
- Auxiliary supply valve flow rates:
 - Standard - 0-25 LPM
 - Optional - 0-13, 0-6 LPM
- Includes 3 main components:
 - A male adapter prob connected to the oxygen outlet.
 - On/Off switch with operating stainless steel lever arm and angled barbed outlet nipple Ø1/4".
 - A Silbermann pressure compensated flowmeter 0-15 LPM or 0-7 LPM.

Silbermann's Quality Management System:

- EN ISO 13485-2016
- ISO 9001-2015
- CE marked





Technical Data

- Operating pressure: 350-500Kpa.
- Pressure compensated 0-15 LPM. Flowmeter calibrated @ 50 psi.
- Flowmeter Material:
 - Body: anodized aluminum*
 - Cover: Transparent non breakable polycarbonate
 - Measuring tube: Transparent non-breakable polycarbonate
 - Valve: brass
 - Orings: Oxygen compatible material
 - Knob: ABS plastics
- Switch-off Valve Material:
 - Body: anodized aluminum
 - Switch: brass
 - Lever Arm: stainless steel

* Brass available on special request



Catalogue Numbers

Type of Adapter*	Silbermann/ Diamond	ENV 737-6	DIN13260-2	Puritan Bennett Type
Wall Mounted	4116100	4116416	4116400	4116150
Rail Mounted**	4116200	4116201	4116202	4116300

* Other adapters on request

** Including: 1m long hose and a universal clamp

Oxygen Switch-off Valve with Auxiliary Instant Supply Valve

Oxygen & Suction
Therapy



Related Products



Humidifier*
P.N - 4101000



Complete DISS
metal nut assembly
for Oxygen
P.N - 2229000



Oxygen DISS,
plastic nut & gland,
for one patient use
P.N - 2261000

* Silbermann's reusable Humidifiers are made of polycarbonate tops, polysulfone jars and are autoclavable at 270°F (132°C).
Jars capacity - 500cc.
Includes safety valve 2 PSI
Inlet - DISS nut assembly.
Includes two flow outlets: 3 - 6mm bushed nipple & 22mm corrugated tube.
ISO color available on request



Assembly - Oxygen Unit



Rail Mounted Oxygen Unit



Oxygen Unit

Distributed By:

SILBERMANN Technologies Ltd.

Address: #5 Harakevet st. POB 4605, Petah-Tikva, 4900715, ISRAEL

Tel: +972-3-9309186 | Fax: +972-3-9315721

E-Mail: info@silbermnn.com, www.silbermnn.com

 **Silbermann**
Medical Gas Systems Since 1950