



We design and manufacture our Terminal Units to last and to ensure a constant and reliable medical gas flow for the patient and the medical equipment.

All outlets are 100 % tested for flow, leakage and are cleaned for oxygen services prior to shipping.

Our outlets are designed by industrial designers to ensure an appealing and pleasant environment for the patient and the medical staff.





Features

- Accepts BS 5682 gas specific adapters.
- All outlets are 100% hydrostatic tested for leaks. The outlets are factory assembled, degreased and cleaned for medical gas service.
- Indexed to eliminate interchangeability of gas services.
- 360° swivel inlet pipe for ease of installation.

Medical Gas System Standards:

 NFPA99, EN ISO 7396-1 & HTM 02-01 requirements.

Applicable Standards:

ISO 9170-1, BS 5682.

Silbermann's Quality Management System:

■ CE marked by the notified body TÜV-SÜD.








Technical Data


- The medical gas outlet shall be SILBERMANN BS type BS 5682
- The outlet shall be gas specific by indexing each front service panel to a gas specific pin indexed rough-in module.
- The outlet of positive and negative pressure shall consist of a primary and secondary check valves to allow testing and on-line maintenance.
- Outlet shall be manufactured with 175 mm long, type K, 12.7 mm (1/2") OD copper pipe. Rough-in pipe shall be silver brazed to the outlet body

- With the rough-in mounted, the outlet shall be depth adjustable to compensate for various wall thicknesses.
- The outlet shall comply with ISO 9170-1 standard typ, quick connect and accept only corresponding BS 5682 gas specific adapter.
- Pressure drop of less than 20 Kpa through the outlet at 120 LPM and 350 Kpa (3.5 Bar) inlet pressure.
- When fully open, the outlet allows a flow of up to 600LPM and pressure of 400 KPa (4 Bars).




Dimensions [mm]

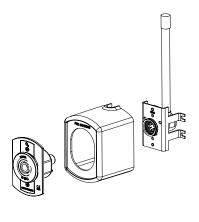

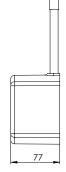



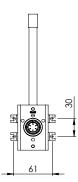



Catalogue Numbers

Gas Service	Color Code	Catalogue No.		
Oxygen	White	3776235		
Medical Air	Black & White	3776257		
Nitrous-Oxide	Blue	3776317		
Vacuum	Yellow	3776318		
Air 800	Black & White	3776319		

DIST	ibute	a B	У.				

SILBERMANN Technologies Ltd.

Address: #5 Harakevet st. POB 4605, Petah-Tikva, 4900715, ISRAEL

Tel: +972-3-9309186 | Fax: +972-3-9315721 E-Mail: info@silbermnn.com, www.silbermnn.com

